
Bed Lever MF 3/11S-CC



Loop with divider rail Handle / Cord and Cleat - resists head traps

Placed between the mattress and the base of your bed. The Bed Lever provides a useful handhold to assist you on and off your bed safely. Correct fitting is essential to the safe use of this product. Take care to read the Use & Fitting Instructions. Also ensure that the Bed Lever is refitted the same way every time the bed is made.

Max user weight	100kg when fitted to a double bed or larger
Loop size	550 H x 225W(81mm between the vertical tubes)
Base size	340 W x 800 L x16 D (all measurements in mm)
Suits bed types	Divan & Inner-sprung beds (all sizes). Does not suit profile beds.
Cleaning	Clean with a soft cloth moistened with soapy water.

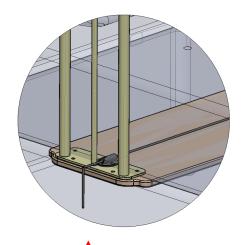




Cord & Cleat can be released by pulling on the cord where it extends under the handle. This means the Bed Lever can be:

1/. Easily removed in an emergency.

2/. Easily refitted each time the bed is made.

Keep the Loop handle away

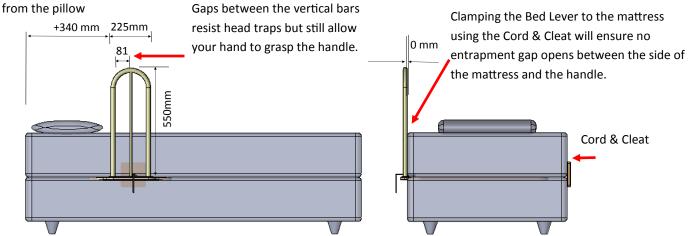

Ph: 0800 102 803, E: reception@multifit.co.nz, Web: www.multifit.co.nz

Multifit Hospital Supplies Ltd 3 bar Bed Lever©270319

Cord & Cleat

9